

TECHNICAL DATA SHEET

Doc No – Ab/TDS/ GLP-169
Rev No: 04 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Powder Funnel, PP

Valid for Catalog No. P70318

SRN (Single Registration Number): IN-MF-000022272

1. Description/Specification

1.1 Description: Powder Funnel , PP

1.1.1 These Powder Funnels are used for transferring material in narrow bottles without spillage.

1.1.2 Various sizes for different container sizes and volumes.

1.1.3 Moulded from high quality virgin Polypropylene to ensure excellent results.

1.2 Design: Powder Funnel , PP

1.3 Performance Specifications:

Code	Internal Dia (mm)	Stem Length (mm)	No. per Case
P70318	65	50	36

1.4 Material:

1.4.1 Polypropylene (PP), Hydrophobic, confirm US FDA 21CFR

1.4.2 Free from Heavy metals, it is safe to use in contact with food stuff, Pharmaceuticals & drinking water.

1.4.3 Product is manufactured without the use of Raw Materials of animal origin. Therefore BSE/TSE issue does not concern this product to the best of our knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

TECHNICAL DATA SHEET

Doc No – Ab/TDS/ GLP-169
Rev No: 04 Date- 05/2024-25

Abdos Labtech Pvt. Ltd. Powder Funnel, PP

1.5 Color: Translucent

1.6 Quality control: Products are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release amendments.

1.7 Basic Features

1.7.1 Made up of High pure virgin Polypropylene (PP)

1.7.2 Temperature range: -80⁰C to +121⁰C

1.7.3 Autoclavability: Yes, Autoclavable at 15 psi and 121⁰C for 20 minutes

1.7.4 Chemical resistance: See website (www.abdoslifesciences.com)

1.7.5 Manufactured in Aseptic Environment

1.7.6 Intended use: Widely used in scientific studies of molecular biology, clinical chemistry, biochemistry and any type of scientific laboratory.

Autoclaving Guideline for Plastics

All plastic consumables are autoclavable at 121⁰C for 20 minutes, 15 psi/1atm. Do not use a dry cycle as this may cause deformation or weakening of the plastic. Allow temperature within the autoclave to return to at least 50⁰C before removing product.

Instruction for Use

Product can be used for food contact also.

Storage temp. Max 30°C

Before using the product check crackness etc.

Check all the information printed on pouches or boxes.

Always wear gloves while handling specimens.

TECHNICAL DATA SHEET

Doc No – Ab/TDS/ GLP-169
Rev No: 04 Date- 05/2024-25

Abdos Labtech Pvt. Ltd.
Powder Funnel, PP


WARNING

DO NOT store strong oxidizing agents.

DO NOT place any plastic labware in a flame.

DO NOT mix any chemicals that may result in a thermic reaction, which can cause product failure

Consult local fire codes prior to storage of flammable liquids.

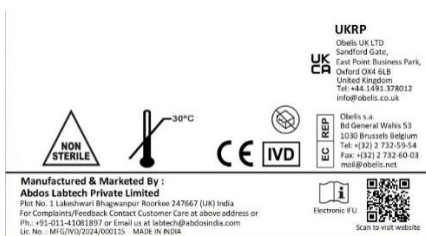

Quality Manager

Abdos Labtech Private Limited

Employee Code 30110

