

“**NEW** ADDITION”

ABDOS
LIFE SCIENCES
Advancing Science with Quality

CompactBLOT Mini Wet Blot System

Abdos introduces Mini Moulded Wet Blotting Transfer System, **CompactBLOT**, an Innovative and Feature rich modular design to run western blotting applications.



Modular Design



Smooth Cassette Locking



Holding Clips with every
Cassette




Single Moulded Tank with
Safety lid & Thumb Locators

COMPACTBLOT - FEATURES

Modular Design :

- Western Blot inner module is designed to be used with the same tank of CompactPAGE vertical electrophoresis.
- The same tank can be utilized for Vertical electrophoresis as well as for western blot transfer applications.

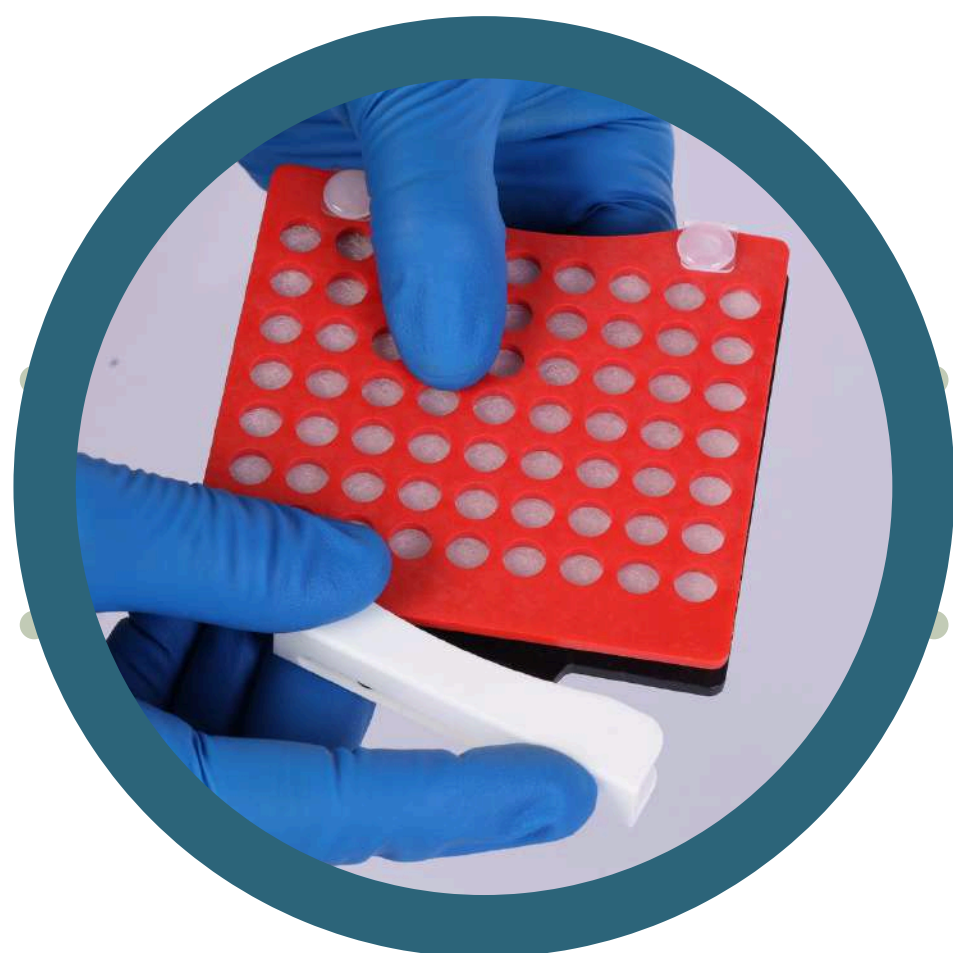


Smooth Safety Lock :

- Each Cassette pair comes with a smooth safety lock to easily open and close the cassette without disturbing the gel and membrane.

Holding Clips :

- The holding clip in the cassette helps the user to easily place or take out the cassettes from the module.

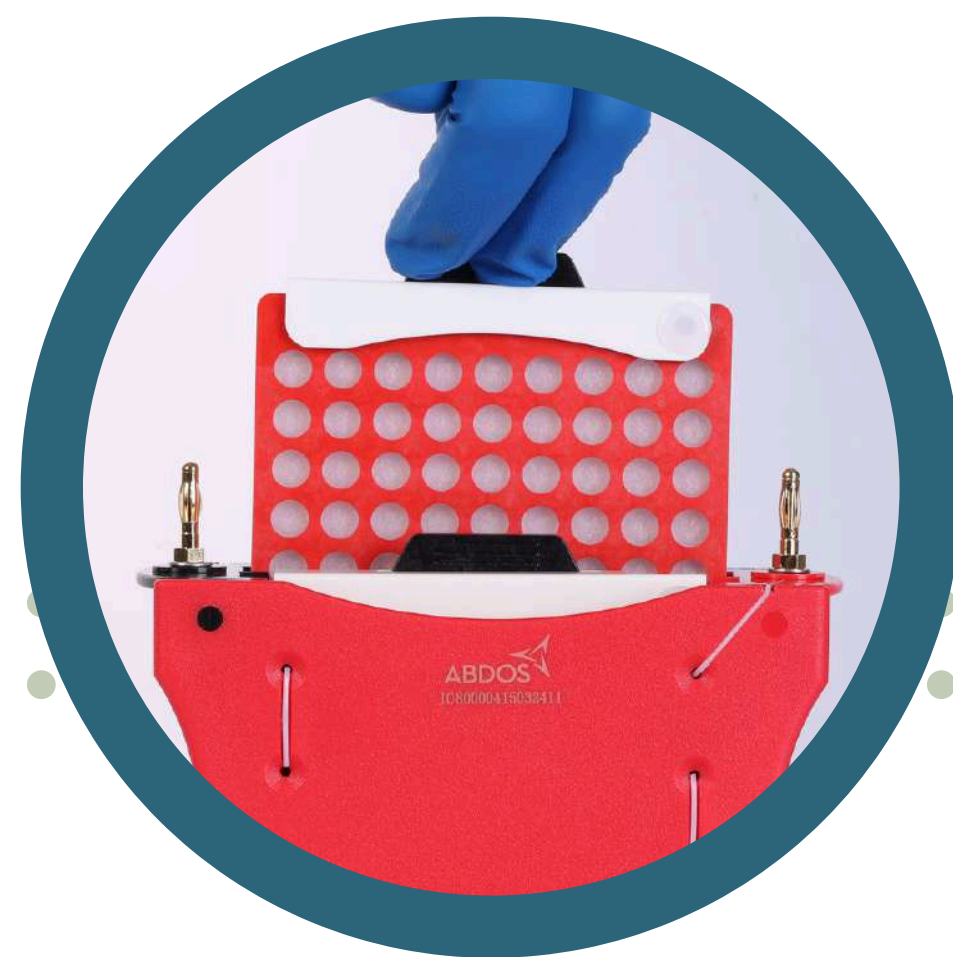

Lid :

- Single moulded lid for long life.
- Easily screw-able cable for safety and storage.
- Vented lid to maintain some amount of heat inside the tank and ensure no condensation on the lid.
- Lid fits only when the orientation is correct.

COMPACTBLOT - FEATURES

Tank :

- Single moulded tank is leakproof, light weight and long lasting.
- Small foot print for easy handling.
- Low buffer requirement design.



Tank base locking :

- Easy locking tank base design to fix the running module so it does not wobble.
- The locking also creates the handling easy for better orientation.

THE COMPACT FAMILY



CompactBLOT
Transfer Module



CompactBLOT
Mini Wet Blot
Complete
System



CompactPAGE
Running Module
& Casting Unit

E12325

COMPACTBLOT MINI WESTERN BLOT

UNIT DIMENSIONS (DXWXH)	10.8x12.5x18.5cm
MATERIAL OF CONSTRUCTION	Polycarbonate
GEL SIZE (WXH)	9 x 7cm
GEL CAPACITY	2 Gels
BUFFER REQUIRED	Maximum - 500ml
MAXIMUM VOLTAGE	200V
RECOMMENDED POWER SUPPLY	E12377

ABDOS CompactBLOT Mini Wet Blot Transfer System - Ordering Information		
Code	Product	No./Case
E12325	CompactBLOT Mini Wet Blot Complete System with CompactBLOT transfer module with cassettes, Running tank, Lid, a pair of Cables and Blotting Papers	1
Spare Accessories		
E12311	CompactBLOT Transfer Module with Cassettes	1
E12310	Wet Blot Transfer Module without Cassettes	1
E12309	Wet Blot Cassette	1